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Qualification of Enterprise

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1. Business Licenses and Certifications

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
Zhejiang Anji Huabu Industrial Co., Ltd.

Shanghai Huabu Medical Technology Co., Ltd.

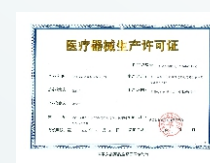Production License of Medical Equipment in Production Base

Production License of Disinfection Product

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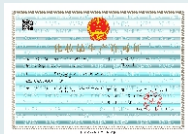
—— (other specifications in accordance with corporate standards or customer orders)

Product Description

Supporting tool for orthopedic wound surgery

Application Scope

Used for fixation of fracture.

Application Method

1. Dipping: At an inclined angle of 45 degrees, dip one roll of products into the water with room temperature, to sip water, until there are no continuous bubbles;
2. Squeezing: Take out the product, and use both hands to squeeze the product from both ends to the middle;
3. Wrapping: Evenly wrap the product around the affected area.
4. Flattening: Flattening by hand at the time of enswathing when enswathing is carried out.

Attentions:


1. This product need to be dipping in the water of 20℃-30℃ about to 5-10 seconds.
2. The wrapping process on one roll of products shall be completed before the curing
3. This product is not in direct contact with the wound.
4. This product is disposable, and it shall be destroyed after being used.
5. Valid Period: 2 years

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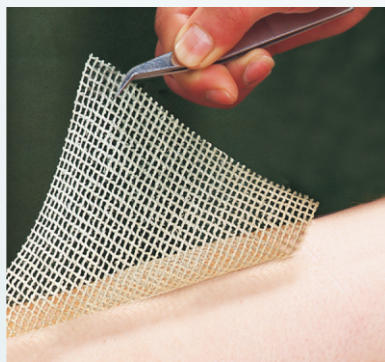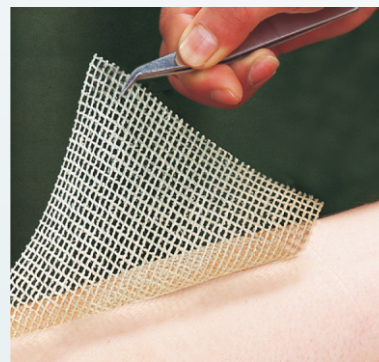
Vaseline Gauze

Specifications

Bagged - flake-shaped product specification

Single-layer, 4-layer, 8-layer, 12-layer, and 16-layer
5×5/7.5×7.5/10×10/10×40/7.5×10/10×20/10×30/15×20/20×20
2.5×90/2.5×20/1.25×180/7.5×45/7.5×90/7.5×22.5/7.5×40/7.5×20
15×90/12.5×22.5/1.2×110-4P/1.25×200/2×200 (cm)

Bagged - rod-shaped product specification

5×200/10×200/15×200/20×200/1×100/2×100/15×100/2.5×460
5×460/1.25×460/0.25×460/1.25×200/2×200/2.5×300/5×300 (cm)

Boxed - flake-shaped product specification:

Single-layer, 4-layer, 8-layer, 12-layer, 16-layer, and 36-layer
5×5/10×10 (cm)

Bagged - folded product specification:

Single-layer, 4-layer, 8-layer, 12-layer, 16-layer, and 36-layer
5×700/10×700/2.5×460/5×460/1.25×460/0.25×460/1.25×200/2×200/
2.5×300 / 5×300 (cm)

—— (other specifications in accordance with corporate standards or customer orders)

Application Scope

This product is applicable to protection, packing and drainage for wound surfaces.

Application Method

1. Please wash your hands before wound disposal.
2. Gently dry the wound and its surrounding skin.
3. Remove the independent packaging bag.
4. Take off the isolation paper on both sides.
Note: Redundant grease is on the isolation paper.
5. Gently put the Vaseline gauze on the wound.
6. If needed, cover the Vaseline gauze swab with the absorbent pad.
7. Fix the gauze with the adhesive tape or bandage.

Attentions

1. In case of serious wound, burn or discomfort arisen from wounds, please immediately go to hospital for treatment.
2. If the packaging bag is damaged, please do not use the product.
3. This product is disposable.
4. In case of wound infection, please immediately go to hospital for treatment.
5. valid Period: 2 years.

This product is based on the actual product. If there is any change, it will not be notified.

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Recommend Products

Zeolite bone Bandage —Patent product

Specifications

7.5*270cm, 10*270cm, 12.5*270cm, 15*270cm, 20*270cm
30*270cm, 7.5*360cm, 10*360cm, 12.5*360cm, 15*360cm
20*360cm, 30*360cm, 7.5*460cm, 10*460cm, 12.5*460cm
15*460cm, 20*460cm, 30*460cm.

—— (other specifications in accordance with corporate standards or customer orders)

Product Description

The zeolite adhesive bandage is made of high-quality pure white natural zeolite powder, calcium sulphate curing agent, non-toxic adhesive material, and high-quality & clean gauze.

Application scope

Used for fixation of fracture. Suitable for patients who need to be orthopedic.

Application Method

1. Dipping: At an inclined angle of 45 degrees, dip one roll of products into the water with room temperature, to sip water, until there are no continuous bubbles.
2. Squeezing: Take out the product, and use both hands to squeeze the product from both ends to the middle, to remove the excess water.
3. Wrapping: Evenly wrap the product around the affected area. Generally, it is suitable to wrap the product for 5-8 layers.
4. Flattening: Within 3 minutes, complete the product shaping and flattening by hand at the time of enswathing, if enswathing is carried out during surgery.

Attentions

1. This product need to be dipping in the water of 20℃-30℃ about to 5-10 seconds.
2. The wrapping process on one roll of products shall be completed before the curing
3. This product is not in direct contact with the wound.
4. This product is disposable, and it shall be destroyed after being used.
5. Valid Period: 24 months
6. Patent NO: ZL20111003.1387.3

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—— (other specifications in accordance with corporate standards or customer orders.)

Product Description

Strip-type or tubular-type, made of elastic or non-elastic materials, not in direct contact with the wound surface.

Application Scope

Used to provide the binding force for wound dressings or limbs and trunks, so as to achieve the purpose of enswathing and fixation.

Application Method

Unpack and take out the product, wrap it around the affected area. The degree of tightness and length shall make the patient feel comfortable.

Attentions

1. This product is not sterilized.
2. This product is disposable, and it shall be destroyed after being used.
3. Please do not use the moldy product.

This product is based on the actual product. If there is any change, it will not be notified.

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5x400cm, 7.5x400cm, 10x400cm, 15x400cm,
2.5x450cm, 5x450cm, 7.5x450cm, 10x450cm,
15x450cm, 20x450cm, 5x270cm, 7.5x270cm,
10x270cm, 15x270cm, 5x500cm, 7.5x500cm,
10x500cm, 15x500cm, 20x500cm

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—— (other specifications in accordance with corporate standards or customer orders)

Product Description

The base fabric of the bandage is made of a combination of a nonwoven fabric and a stretch yarn material. The adhesive layer of the bandage is composed of a pressure-sensitive adhesive or a hot-melt adhesive non-toxic sanitary material. This product contains elastic or non-elastic materials, and it should not in direct contact with the wound surface.

Application Scope

Used to provide the binding force for wound dressings or limbs and trunks, so as to achieve the purpose of enswathing and fixation.

Application Method

1. Unpacking: due to the feature of product: it should not be damaged or leaked before open it.
2. Enswathing: Used for protection for elbows:
 - (1) gently make the arm close to the shoulder and form an angle of 90° paste the bandage along the elbow from the inner side of the wrist to the shoulder.
 - (2) Bend the elbow, lift the wrist, and start to paste the bandage from the inner side of the upper wrist; stretch out the elbow, and then paste the bandage on the skeletal muscle of the wrist.
3. Used for stretch of wrist acupoint:
 - (1) Cut the adhesive tape into segments and embed them into the gap between fingers.
 - (2) One is attached to the wrist from the back of the hand.
 - (3) The other is attached to the wrist from the outstretched palm.
 - (4) In order to avoid that the wrist is oppressed, please wind once again.

This product is based on the actual product, if there is any change, it will not be notified.

Attentions

This product is not sterilized; This product is disposable, and it shall be destroyed after being used.

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Recommend Products

Cast Padding

(Whole cotton / chemical fiber / bamboo fiber)

Specification

5 × 270cm, 7.5 × 270cm, 10 × 270cm, 15 × 270cm,
10 × 300cm, 15 × 300cm, 10 × 450cm, 20 × 270cm,
5 × 360cm, 7.5 × 360cm, 10 × 360cm, 15 × 360cm,
20 × 360cm, 4 × 300cm, 6 × 300cm, 8 × 300cm,
20 × 300cm, 10 × 200cm, 5 × 450cm, 7.5 × 450cm,
15 × 450cm, 20 × 450cm

—— (other specifications in accordance with corporate standards or customer orders)

Application Scope

General protection for patients, so as to avoid the injury arisen from other devices or external factors.

Application Method

1. Unpack the self-sealing bag, and take out the cast padding.
2. Evenly wrap it around the affected area.
3. Pull apart the padding, and keep the surface flat.

Attentions

1. When enwinding the cast padding, please do not use too much force, so as to avoid tearing or breakage.
2. This product is not sterilized.
3. This product is disposable, and it shall be destroyed after being used.

Production Description


General protective article or material for patients during treatment.

This product is based on the actual product. If there is any change, it will not be notified.

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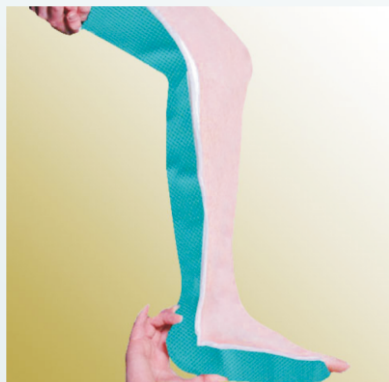7.5x30cm, 7.5x90cm, 10x40cm, 10x50cm,
10x75cm, 12.5x76cm, 12.5x115cm, 15x75cm,
15x115cm, 10x38cm, 10x76cm, 12.5x76cm,
12.5x114cm, 7.5x88cm, 12.5x114cm, 15x76cm,
15x114cm

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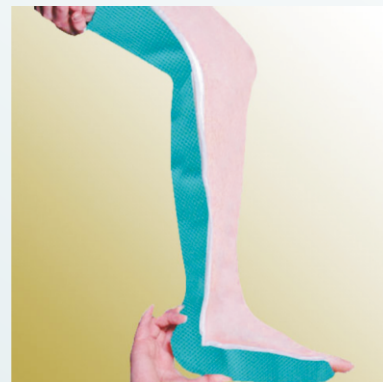
—— (other specifications in accordance with corporate standards or customer orders)

Production Description

Supporting tool for orthopedic trauma surgery.

Application Scope

Used for external fixation of fracture.

Application Method

1. The operator should wear gloves before the operation to avoid contamination of hands by glue.
2. Select appropriate specification and quantity of the product based on the actual needs of the wound of the patient, and unpack the product one by one during operation.
3. According to use needs, before use, take out the unpacked product, hold the splint by hand, and then immerse it into the water with temperature of 20-45℃ (68-75°F) for 3-10 seconds until the product is fully soaked.
4. Take out the product from the water, hold the splint with both hands, slightly use more force, and squeeze out the excess water in the product.
5. Put the soaked and dried product on the patient's wound; shape the product within 3-8 minutes.
6. Fetch one roll or several rolls of gauze bandage or cotton pad; enwind them onto the splint in a certain sequence and direction, in order to fix the

This product is based on the actual product. If there is any change, it will not be notified.

Attentions

1. As the product is subject to gradually curing in contact with air, please do not use the product that remains unused after unpacking, has leaks in package, has damaged package due to inappropriate storage or is expired in shelf-life.
2. Operators must put on protective gloves before the use of the product to avoid adhesion of the glue to hands.
3. This product is not in direct contact with the wound.
4. This product is disposable.
5. Valid period: 2 years

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—— (other specifications in accordance with corporate standards or customer orders)

Product Description

Strip-type or tubular-type, made of elastic or non-elastic materials, not in direct contact with the wound surface.

Application Scope

Used to provide the binding force for wound dressings or limbs and trunks, so as to achieve the purpose of enswathing and fixation.

Application Method

1. The operator should wear gloves before the operation to avoid contamination of hands by glue.
2. Select appropriate specification and quantity of the product based on the actual needs of the wound of the patient, and unpack the product one by one during operation.
3. According to use needs, before use, take out the unpacked product, hold the splint by hand, and then immerse it into the water with temperature of 20-45℃ (68-75°F) for 3-10 seconds until the product is fully soaked.
4. Take out the product from the water, hold the splint with both hands, slightly use more force, and squeeze out the excess water in the product.
5. Put the soaked and dried product on the patient's wound; shape the product within 3-8 minutes.
6. Fetch one roll or several rolls of gauze bandage or cotton pad; enwind them onto the splint in a certain sequence and direction, in order to fix the position of the splint.

Attentions

1. As the product is subject to gradually curing in contact with air, please do not use the product that remains unused after unpacking, has leaks in package, has damaged package due to inappropriate storage or is expired in shelf-life.
2. Operators must put on protective gloves before the use of the product to avoid adhesion of the glue to hands.
3. This product is not in direct contact with the wound.
4. This product is disposable.
5. Valid period: 2 years

This product is based on the actual product. If there is any change, it will not be notified.

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—— (other specifications in accordance with corporate standards or customer orders)

Application Scope

Used to provide the binding force for wound dressings or limbs and trunks, so as to achieve the purpose of enswathing and fixation.

Application Method

1. Unpack, and take out the product.
2. Wrap it around the affected area, The degree of tightness and length shall make the patient feel comfortable.

Attentions

1. This product is not sterilized;
2. This product is disposable, and it shall be destroyed after being used.
3. Valid Period: 2 years

Production Description

Strip-type or tubular-type, made of elastic or non-elastic materials, not in direct contact with the wound surface.

This product is based on the actual product. If there is any change, it will not be notified.

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1.5 × 200cm/卷, 2 × 200cm/卷
2.5 × 200cm/卷, 1.5 × 100cm/卷
2 × 100cm/卷, 2.5 × 100cm/卷

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